

GOLDEN
Passport[™]
HI LOW ADJUSTABLE BED

Have the *Rest* of Your Life

IN A GOLDEN PASSPORT[™] HI LOW ADJUSTABLE BED

 **GOLDEN**
POSITIONING YOU FOR LIFE



Rich Golden pictured with his children, third generation Golden Team Members! Pictured left to right: Rich Jr., Alyssa, Rich & Mike.

A MESSAGE FROM RICH GOLDEN, CEO

Golden is known for comfort, quality and innovation. Founded in 1985, we manufacture power lift recliners, scooters and power wheelchairs that push the boundaries of expectation. We are proud to re-enter the bedding space by adding the new **Golden Passport™ Hi Low Adjustable Bed** to our product line.

The Golden Passport will provide users with peace of mind knowing they can pursue their own independence. Whether customers are seeking the many health benefits of a hi low adjustable bed, recovering from an injury or simply have mobility concerns, the Golden Passport gives users (and caregivers) an added sense of security and peace of mind in their home.

OUR GOAL is to provide YOU with the best "REST" of your life.

A handwritten signature in black ink that reads "Rich Golden".

OUR VISION...

To build the most innovative and comfortable products that put you and your family's lives in motion.



GOLDEN
Passport™
HI LOW ADJUSTABLE BED

Understanding the Wellness Benefits of an Adjustable Bed

Did you know that the way you sleep, and the amount of support received during your time of rest makes a big difference in your overall health?

When you have the option to change positions while resting, you can combat many different health problems, specifically – pains and ailments.

With an adjustable bed, it is possible to alter your sleep position precisely as you wish by independently inclining or declining your head and/or legs into a position that more closely matches the contours of your body. The more points of adjustments you have, the greater the benefits to your health. You can achieve many wellness positions like Zero-Gravity and the Trendelenburg position.

Benefits can include:

- **Improved sciatica** by supporting the base of the spine
- **Relief of pressure** and compression in the back
- **Reduction in snoring** & improved quality of sleep
- **Lessened symptoms** of some types of **sleep apnea**
- **Healthier digestion** with a 6" recommended incline
- **Optimized circulation & reduced swelling** by elevating the lower limbs

What is a "Hi Low" Adjustable Bed?

A "Hi Low" bed is height adjustable via remote control, with ability to be custom elevated or lowered, based on the user's or caregiver's specific needs.



The Golden Passport™

HI LOW ADJUSTABLE BED WITH DUAL ZONE MASSAGE



The Golden Passport offers many unique and luxurious properties, and is available in Twin XL, Queen, and Split King (Dual Twin XL) & one standard foam mattress.


Unique Mattress Design for Optimal Flexibility

The Golden Passport's mattress provides comfort and support to ensure a restful night's sleep. Our memory foam mattress is made of three layers of foam for solid support paired with the soft luxurious feel of memory foam. It features a **unique transfer edge**, a 3" firm block of foam that runs the length of the mattress, which provides stability for safe transfers when entering and exiting the bed.

The 2" top layer is GelFlex foam that helps to keep you cool while you sleep.

The 2" middle layer is firmer than the top layer with great rebound to enhance the supportive feel of the mattress.

The 6" base layer is made of a very firm core foam that is engineered to contour perfect to the Golden Passport's adjustable base, making sleep more comfortable than other memory foam mattress. You'll enjoy firm, full body support that allows you to easily turn over at night and helps to maintain your mobility when transferring in and out of bed.

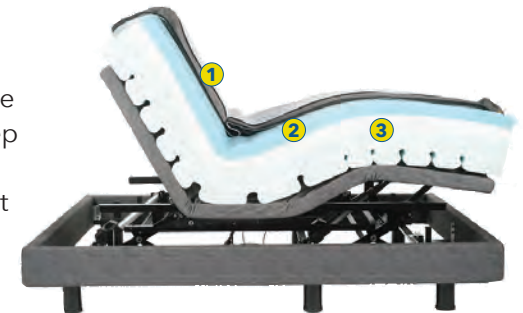

HI LOW HEIGHT OPTIONS



LEG EXTENSIONS FOR THE MOST HEIGHT FLEXIBILITY



TRANSFER EDGE



NUMEROUS POSITIONING OPTIONS



WIRELESS REMOTE



LOCKING FOOT CASTERS



SIDE RAILS WITH & WITHOUT COVERS



Bed Frame Features:

- 880 lb. lift capacity
- 3-piece leg extensions included (2", 2.5" & 5") for varied height options
- Up to 11" vertical rise
- **Lowest Position:** 10.5" from floor to top of deck (*no leg sections added*)
- **Highest Position:** 28" from floor to top of deck (*with all leg sections added*)
- Dual-Zone Vibrating Massage with timer
- Feature-rich, 22-button wireless remote with two programmable buttons and four pre-programmed positioning buttons including TV, Zero Gravity, Flat & Quiet Sleep (a position that helps open the airways by slightly elevating the head)
- Tilt feature in the head section to achieve reverse Trendelenburg position
- Wall hugger design
- Convenient lighting features include underbed LED lighting and flashlight in remote
- Headboard bracket
- Front and side retainers to hold the mattress securely in place
- Battery backup
- Split head & foot bases for easier & lighter weight delivery

Bed Frame & Mattress Sizes Available:

- Twin XL (38" x 80")
- Queen (59" x 80")
- Split King (76" x 80" - 2 x Twin XL)

Limited Lifetime Warranty:

- Lifetime on Frame
- 1 Year on Mechanical & Electrical Labor
- 3 Years on Mechanical & Electrical Parts
- 5 Years on Mattress

Accessories Sold Separately:

- Side Rails (Set of 2)
- Side Rail Covers (Set of 2)
- Locking 2.5" Foot Casters (Set of 8)

Mattress Features:

- 350 lb. weight capacity
- High density base foam for maximum durability and resiliency
- SmartFoam for optimal head & foot articulation
- Cooling gel memory foam
- Premium-grade, machine washable cover with stretch fibers for a supremely comfortable feel
- Firm transfer edge (or edge guard) to stay centered while resting and assist with safe transfers in and out of bed
- Premium packaging includes shrink wrapped mattress in a manageable size

